



Discover your protein signature

State-of-the-art benchtop solution for multiplex protein detection



Olink® Signature Q100

- Dedicated benchtop solution designed for readout of Olink's high quality Target 96, Target 48, Flex and Focus biomarker panels
- User-focused design and intuitive interface, including integrated software and IFC loader
- Small, compact footprint
23.6 in (60 cm) x 10.6 in (27 cm) x 21.6 in (55 cm)
- Low investment threshold
- NPX Signature software for intuitive data QC and normalization.

www.olink.com

© 2025 Olink Proteomics AB, part of Thermo Fisher Scientific.
Olink products and services are For Research Use Only. Not for use in diagnostic procedures.
This document is not intended to convey any warranties or representations of any kind.
1320, v1.2, 2025-10-09

Olink Signature Q100

Supported product lines

Olink® Target 96

Includes 13 carefully designed panels for human (12) and mouse (1) representing a library of over 1100 proteins. Each panel is based on a specific disease area or biological process, and enables the relative quantification of 92 proteins across 88 samples simultaneously, using just 1 µL per sample.

Olink® Target 48

Delivers relative or absolute quantification of up to 45 carefully selected proteins across 40 samples simultaneously from just 1 µL of sample. Panels available for human (3) or mouse (1). Panels can be combined and processed simultaneously for a comprehensive analysis.

Olink® Flex

Build your own biomarker panel. Freely mix and-match 5 - 30 proteins from our library of around 200 pre-validated inflammatory focused proteins. Delivers relative or absolute quantification.


Olink® Focus

A custom mid-plex solution targeting 21 customer selected proteins, with a choice of absolute or relative quantification.



Olink® Target Panels

Olink® Target 96 Cardiometabolic	Olink® Target 96 Inflammation	Olink® Target 96 CVD II & III	Olink® Target 96 Metabolism
Olink® Target 96 Development	Olink® Target 96 Immune Response	Olink® Target 96 Immuno-Oncology	Olink® Target 96 Oncology II & III
Olink® Target 96 Mouse Exploratory	Olink® Target 96 Neurology	Olink® Target 96 Neuro Exploratory	Olink® Target 48 Immune Surveillance
Olink® Target 48 Neurodegeneration	NEW Olink® Target 48 Cytokine	Olink® Target 48 Mouse Cytokine	

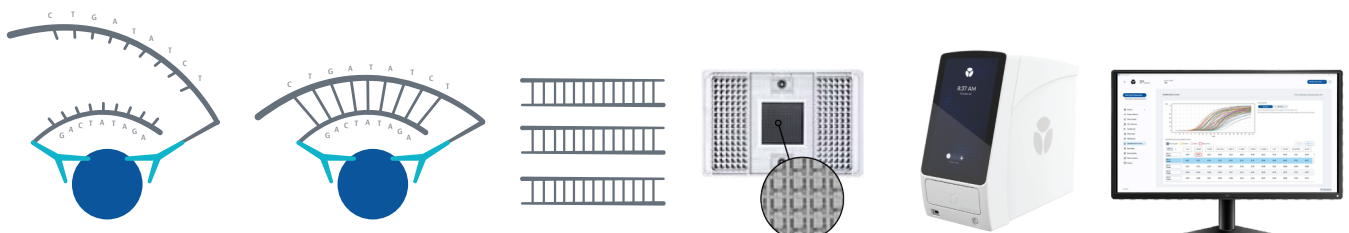

Lab workflow for Olink Signature Q100. Antibody pairs labeled with DNA oligonucleotides bind target antigen in solution, allowing hybridization and extension by DNA polymerase. This newly created piece of DNA barcode is amplified by standard PCR and then transferred to an integrated microfluidic chip (IFC), which is loaded into the Olink Signature Q100 instrument for qPCR and data readout.